


OLYMPUS®

Olympus is about life. About photographic innovations that capture precious moments of life. About advanced medical technology that saves lives. About information- and industry-related products that make possible a better living. About adding to the richness and quality of life for everyone. Olympus. Quality products with a **FOCUS ON LIFE**

HOFFMAN MODULATION
CONTRAST UNIT

IX-HMC

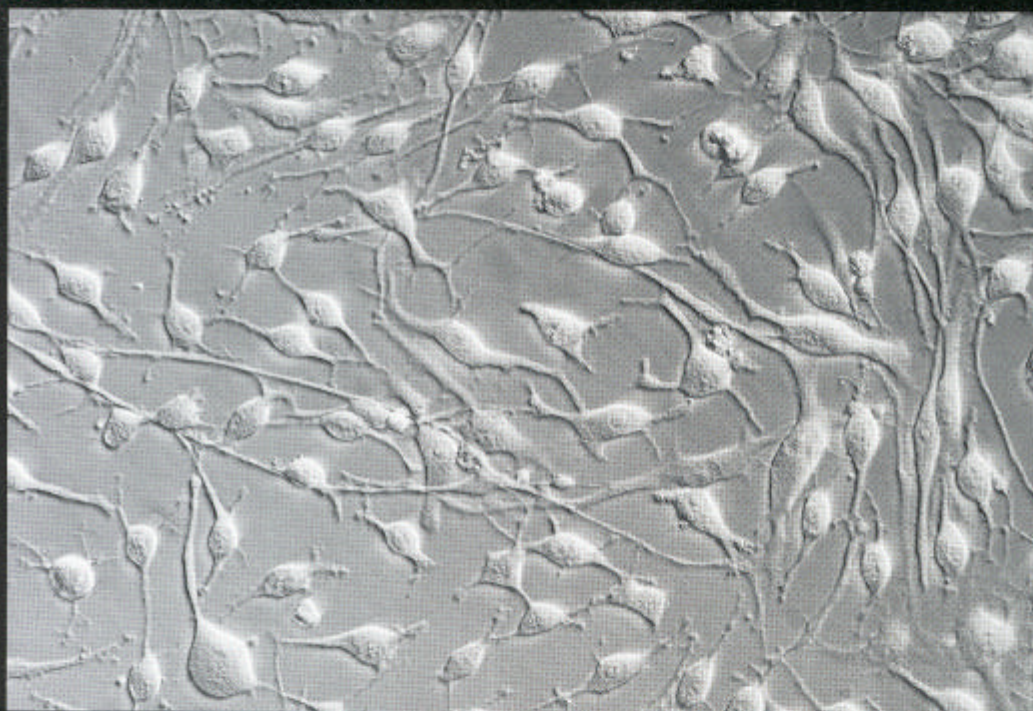

Hoffman Modulation
Contrast®



Nomarski DIC

The image at right was produced
using a plastic chamber

Plastic chamber eases observation of 3D-like contrast.

Hoffman Modulation Contrast[®] unit


- Objectives, condensers (W.D. 40mm, N.A. 0.5, 4-holder) and polarizers developed specifically for Hoffman Modulation Contrast.
- Simple, clear 3D-like contrast images for easier manipulation.
- Plastic chambers, not applicable with Nomarski DIC, are possible to use.
- Thick specimens which cannot be observed with phase-contrast due to their halos can be observed with Hoffman Modulation Contrast.



Hoffman Modulation Contrast objectives



Item	N. A.	W. D.
HMC10 LWD CPLAN	0.25	9.8mm
HMC20 LWD LCA	0.40	3.2mm
HMC40 LWD LCA	0.55	1.2mm

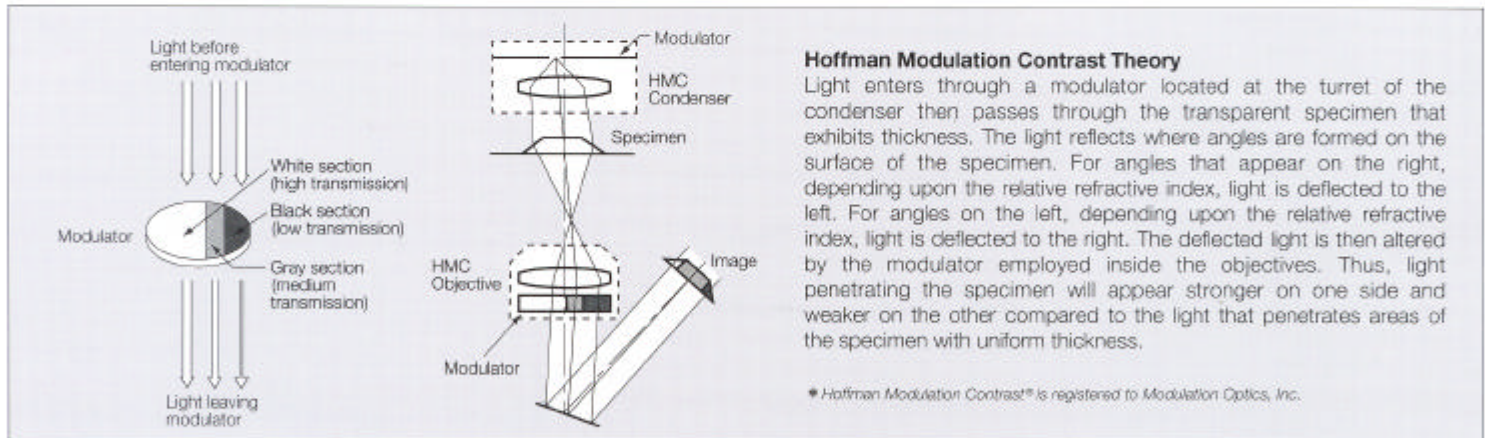

Specifications are subject to change without any obligation on the part of the manufacturer.

Olympus business areas



OLYMPUS[®]

FOCUS ON LIFE

OLYMPUS OPTICAL CO., LTD.
2-43-2 Hatagaya, Shibuya-ku, Tokyo, Japan
OLYMPUS OPTICAL CO. (EUROPA) GMBH
Postfach 10 49 08, 20094, Hamburg, Germany
OLYMPUS OPTICAL CO. (U.K.) LTD.
2-B Honduras Street, London EC1Y 0TX, United Kingdom.